


PARAFFIN OVENS UNpa

according to DIN 12880:2007-05; EN 61010-1 (IEC 61010-1); EN 61010-2-010

Standard units are safety-approved
and bear the test marks:

- Interior: stainless steel, material 1.4301 (ASTM 304), with all-round deep-drawn ribs to integrate the large-area heating with ceramic-metal sheath
- Housing: textured stainless steel 1.4016 (ASTM 430), rear zinc-plated steel 1.0037 (St37-2), intuitively operated TwinDISPLAY (TFT colour display) with touchscreen
- Door: fully insulated stainless steel doors with 2-point locking (compression door lock)
- Connection: mains cable with plug (German type)
- Installation: 4 adjustable appliance feet
- Interfaces:

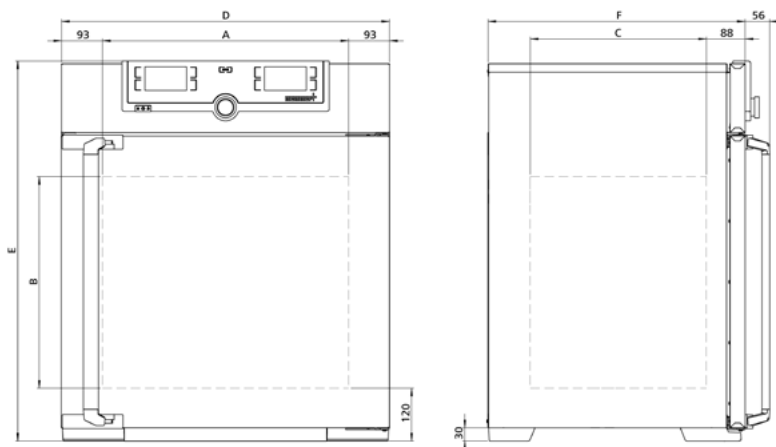

Model sizes/Description				30	55	75	110	160
Stainless steel interior	Volume	l		32	53	74	108	161
	Width	(A) mm			400		560	
	Height	(B) mm		320	400	560	480	720
	Depth	(C) mm		250		330	400	
	Max. number of grids/shelves	number		3	4	6	5	8
	Max. loading per grid/shelf	kg				20		
	Max. loading of chamber	kg		60	80	120	175	210
	Max. loading per slide-in drip tray	kg			1.5		3	
	Max. loading per bottom drip tray	kg			1.5		3	
Textured stainless steel exterior	Width	(D) mm			585		745	
	Height	(E) mm		704	784	944	864	1104
	Depth (without door handle, depth of handle +56 mm)	(F) mm		434		514	584	
Standard equipment	Standard works calibration certificate (measuring point chamber center)	°C				+80		
	Stainless steel grid, electropolished	number			1		2	
Temperature	Working temperature range	°C				at least 5 above ambient temperature up to +80		
	Setting temperature range	°C				+20 to +80		
	Setting accuracy	°C				0.1		
Further data	Electrical load at 230 V, 50/60 Hz	W		1600	1000	1250	1400	1600
	Electrical load at 115 V, 50/60 Hz	W		800		900		
Packing data	Net weight	kg		45	55	66	75	96
	Gross weight (packed in carton)	kg		61	74	85	100	122
	Width	mm		660		730	830	
	Height	mm		890	950	1130	1050	1300
	Depth	mm		650		670	800	
Order No. Paraffin Ovens				UN30pa	UN55pa	UN75pa	UN110pa	UN160pa